

mediTECH

MEDICAL GAS EQUIPMENT

microVENT™ Utilityvent


Emergency Automatic Resuscitator Kit



- Automatic resuscitation.
- Minimal product training needed.
- Easy to use in the dark or with poor access.
- Pre-set for adult casualties for ease of use in Industrial conditions.
- Suitable for use in non-breathable atmosphere.



microVENT™ Utilityvent Kit
Part No. 670-0340-01

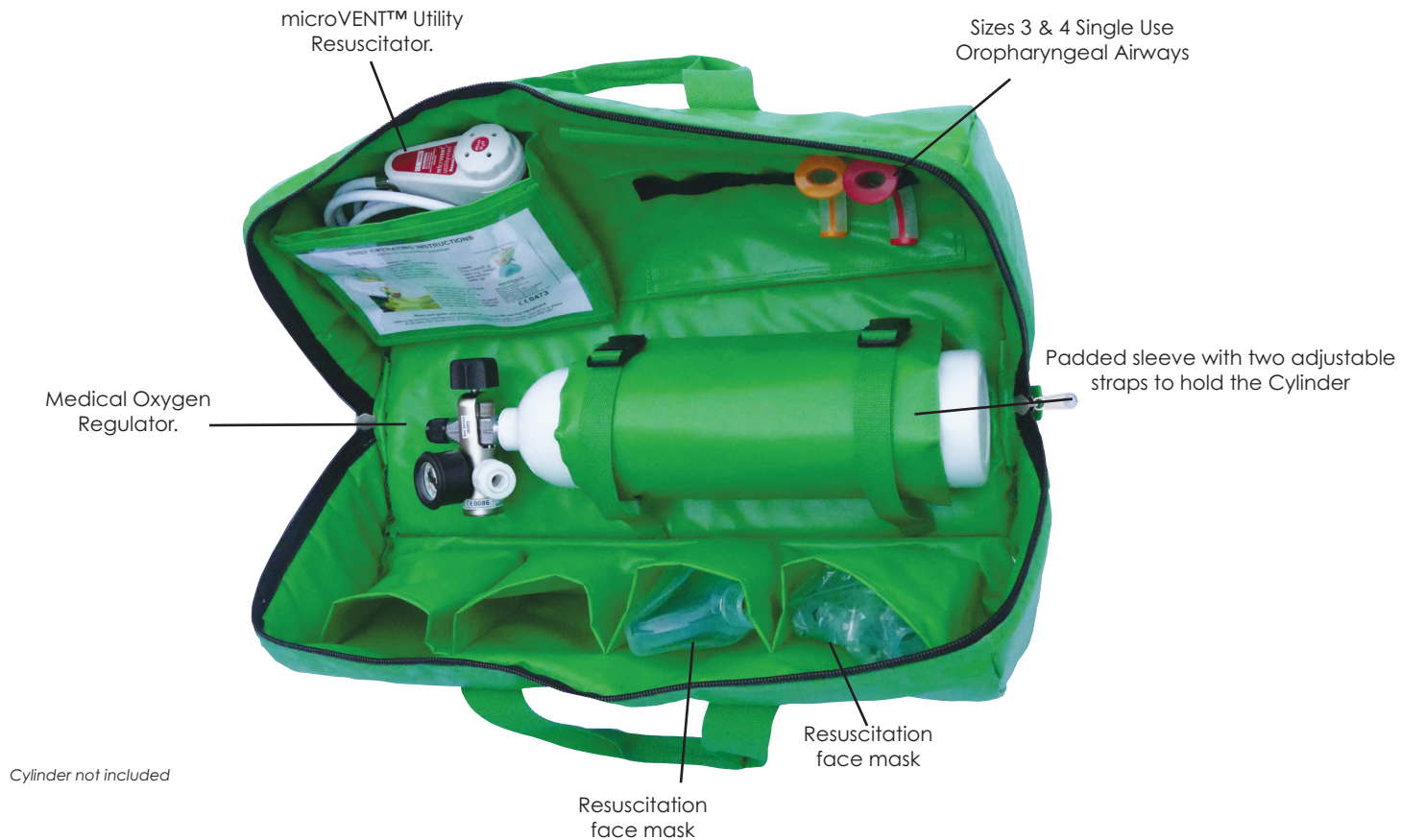

The microVENT Utilityvent kit is designed for use in an emergency when breathing has stopped, or when breathing assistance is needed. At the heart of the kit is the microVENT Utility Resuscitator. This is an oxygen powered resuscitator that is pre-set to breathe for the adult casualty - 600ml breaths, 12 times a minute. When in use, if the casualty inhales, the Utilityvent supplies a synchronised breath of 100% oxygen. The brass safety oxygen regulator fits to pin-index medical oxygen cylinders. The microVENT case is reinforced and padded to hold and protect the equipment. It is fitted with carrying handles and a shoulder strap. The resuscitator can be swiftly carried to the casualty.



microVENT™ Utilityvent Resuscitator

Part No. 670-0339-00

A hand held, light weight Resuscitator with compact dimensions, the microVENT™ Utilityvent has ergonomic styling to give easy operation with controls in a natural position. It is constructed with a tough reinforced plastic casing to withstand use in difficult situations.

The Utilityvent's ergonomic styling gives easy operation of the resuscitator. The patient valve features 22/15 mm connections to take standard airway devices and face masks. The oxygen supply hose is two metres long and made of resilient reinforced plastics to retain flexibility in cold conditions and is permanently attached to the regulator, ready for immediate use.

- Oxygen gauge is luminous to aid use in low light.
- microVENT is used by emergency services around the world.



Medical Oxygen Pin Index Regulator

Part No. 55B 678-0001-10

The microVENT™ Mariner Kit includes two Meditech Oxygen Therapy Regulators with a Pin Index Inlet connector and a BS Medical Oxygen Schrader outlet



medITECH

MEDICAL GAS EQUIPMENT

B.N.O.S. Meditech Ltd.

**9 Fifth Avenue, Bluebridge Industrial Estate,
Halstead, Essex CO9 2SZ, England, UK.**

Tel: +44 (0)1787 479475 • Fax: +44 (0)1787 477747

E-mail: sales@meditech.uk.com • www.meditech.uk.com

TO PURCHASE OR FOR FURTHER DETAILS:

**Contact us by email or telephone for more details,
including specific part numbers, prices, etc.**

B.N.O.S. Meditech Ltd.
9 Fifth Avenue, Bluebridge Industrial Estate,
Halstead, Essex CO9 2SZ, England, UK.
Tel: +44 (0)1787 479475 • Fax: +44 (0)1787 477747
E-mail: sales@meditech.uk.com • www.meditech.uk.com



Components manufactured by BNOS Meditech and supplied within the kit are CE marked in accordance with the European Medical Devices Directive 93/42/EEC and UKCA marked accordance with the UK Medical Devices Regulations 2002. All items in this range are sold with a return to manufacturer warranty against manufacturing defects. These products are part of a series of medical gas and resuscitation products available from B.N.O.S. Meditech. If the product you require is not shown here, or if you require more information, please contact your distributor or contact us directly. Specifications may be subject to change without notice.

B.N.O.S. Meditech Ltd Terms and Conditions of sale apply, available on application or from our website.

Distributor: