



## 4 Way Flow Selector and Frame

### BS Schrader Probe on 2M Supply Hose

- **No Calibration required EVER**
- **Lightweight construction**
- **Available with International Fittings**
- **7 Year Warranty**
- **7 Year service Interval**
- **Mass Casualty Device**

Excellent performance and robustness, coupled with a very lightweight design. Manufactured from a strong aluminium frame with folding carry handle, extending hooks for hanging on the side of a stretcher and plastic feet on the base. Supplied with 2M supply hose to connect to 4 bar oxygen source (longer hose length available upon request)

For use in all Hospital and Emergency Environments, the 4 Way Flow Selector is designed to be plugged into Medical Oxygen Wall Outlets.

The Meditech 4 Way Flow Selector, has been designed to deliver precision flow rates, without the need for any calibration.

Supplied with Push Fit "fir-tree" for direct attachment of Oxygen Therapy accessories.

**AVAILABLE GAS TYPES:**  
OXYGEN | MEDICAL AIR | ENTONOX

#### THERAPY OUTLET TYPES



FIR TREE



9/16<sup>th</sup> THREADED THERAPY CONNECTOR



3/8" BSP THREADED THERAPY CONNECTOR




1/4" BARB

#### PROBE TYPE



BS PROBE

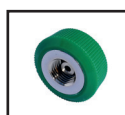

DISS



DIN



AGA



AFNOR



SIS



CZECH



UNIFOR



# SPECIFICATION

Intended Gas: .....	Medical Oxygen (O <sub>2</sub> )
Inlet pressure: .....	4 bar (400 kPa)
Set Flow Rates in lpm: .....	Three flow rates options available: 1/2, 1, 2, 3, 4, 6, 8, 10, 15 lpm 1/2, 1, 2, 3, 4, 6, 8, 10, 15, 25 lpm 10, 20, 30, 35, 40, 50, 55, 60, 70 lpm
Inlet Connections: .....	Schrader (BS 5682:2015) and other international connections
Outlet (Patient) Connections. ....	Firtree style barbed outlet.
Intended Life: .....	10 Years
Service Interval: .....	7 Years
Warranty Period: .....	7 Years
Complies With: .....	ISO 15002: 2008
Weight including hose:.....	3.4Kg
Dimensions: .....	Height: 18.5CM x Length: 46CM Height with handle raised: 24CM Depth: 15.5CM

## TO PURCHASE OR FOR FURTHER DETAILS:

Contact us by email or telephone for more details, including specific part numbers, prices, etc.

**B.N.O.S. Meditech Ltd.**  
**9 Fifth Avenue, Bluebridge Industrial Estate,**  
**Halstead, Essex, CO9 2SZ, England, UK.**  
**TEL: +44 (0) 1787 479 475**  
**FAX: +44 (0) 1787 477747**  
**EMAIL: [sales@meditech.uk.com](mailto:sales@meditech.uk.com)**  
**WEBSITE: [www.meditech.uk.com](http://www.meditech.uk.com)**



**medITECH**  
MEDICAL GAS EQUIPMENT

Distributor:

Meditech Flow Selectors are CE marked in accordance with the European Medical Devices Directive 93/42/EEC and UKCA marked accordance with the UK Medical Devices Regulations 2002. All items in this range are sold with a return to manufacturer warranty against manufacturing defects. These products are part of a series of medical gas and resuscitation products available from B.N.O.S. Meditech. If the product you require is not shown here, or if you require more information, please contact your distributor or contact us directly. Specifications may be subject to change without notice.

B.N.O.S. Meditech Ltd Terms and Conditions of sale apply, available on application or from our website.

Doc 771M Issue D 28/08/2024