

Straight Bar, Rail Mounted Flow Selector 2M Supply Hose, BS Probe -100mm



- **No Calibration required EVER**
- **Lightweight construction**
- **Available with International Fittings**
- **7 Year Warranty**
- **7 Year service Interval**

Excellent performance and robustness, coupled with a very lightweight design.

For use in all Hospital and Emergency environments, the Straight Bar Flow Selector and Outlet is designed to be plugged into Medical Oxygen Wall Outlets, whilst being mounted on a medical rail.

The Meditech Straight Bar Flow Selector and Outlet has been designed to deliver precision flow rates, without the need for any calibration. Its fitted with a universal size rail mount for most rail systems. It has a standard 2M hose length (other hose lengths available upon request).

Can be supplied in either of two formats: Push-Fit 'Firtree' Fitting for direct attachment on a Oxygen Therapy delivery system With a threaded outlet fitting for the attachment of a humidifier.

THERAPY OUTLET TYPES



FIR TREE



9/16th THREADED THERAPY CONNECTOR



3/8" BSP THREADED THERAPY CONNECTOR




1/4" BARB

AVAILABLE GAS TYPES:
OXYGEN | MEDICAL AIR | ENTONOX

PROBE TYPE



BS PROBE

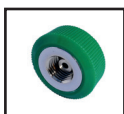

DISS



DIN



AGA



AFNOR



SIS



CZECH



UNIFOR



OUTLET TYPES



BS SCHRADER



DISS



DIN




AGA



AFNOR




SIS

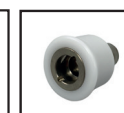

CZECH



UNIFOR

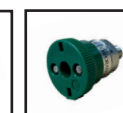

OHIO



SPECIFICATION

Intended Gas:	Medical Oxygen (O ₂ , MA4, O ₂ /N ₂ O)
Inlet pressure:	4 bar (400 kPa)
Set Flow Rates in lpm:	Three flow rates options available: 1/2, 1, 2, 3, 4, 6, 8, 10, 15 lpm 1/2, 1, 2, 3, 4, 6, 8, 10, 15, 25 lpm 10, 20, 30, 35, 40, 50, 55, 60, 70 lpm
Inlet Connections:	Schrader (BS 5682:2015) and other international connections
Outlet (Patient) Connections.	Firtree style barbed outlet, 9/16 threaded outlet, 3/8" BSP threaded outlet & 1/4" barb
Intended Life:	10 Years
Service Interval:	7 Years
Warranty Period:	7 Years
Complies With:	ISO 15002: 2008

TO PURCHASE OR FOR FURTHER DETAILS:
Contact us by email or telephone for more details, including specific part numbers, prices, etc.

B.N.O.S. Meditech Ltd.
9 Fifth Avenue, Bluebridge Industrial Estate,
Halstead, Essex CO9 2SZ, England, UK.
Tel: +44 (0)1787 479475 • Fax: +44 (0)1787 477747
E-mail: sales@meditech.uk.com • www.meditech.uk.com



medITECH
MEDICAL GAS EQUIPMENT

Distributor:

Meditech Flow Selectors are CE marked in accordance with the European Medical Devices Directive 93/42/EEC and UKCA marked in accordance with the UK Medical Devices Regulations 2002. All items in this range are sold with a return to manufacturer warranty against manufacturing defects. These products are part of a series of medical gas and resuscitation products available from B.N.O.S. Meditech. If the product you require is not shown here, or if you require more information, please contact your distributor or contact us directly. Specifications may be subject to change without notice.

B.N.O.S. Meditech Ltd Terms and Conditions of sale apply, available on application or from our website.

Doc 785M Issue D 28/08/2024